

New!

Intensiv Proxoshape Flexible

Clinically tested by the University of Zurich, Switzerland

Oscillating diamond-coated file for the removal of excess material and contouring of restored proximal surfaces.


Benefits

- Simplified access to the interdental space
- Improved adjustment to the proximal filling morphology
- No iatrogenic damage to adjacent tooth due to one side diamond-coating
- Prevention of irregular surface waves created by use of rotating instruments

The Intensiv Proxoshape Flexible file is, due to its minimal file thickness of 0.1 mm, so thin that difficult to access interdental spaces can be treated, Intensiv Proxoshape Flexible allows easier access to the interdental space. In addition, the thin file is sufficiently flexible to adapt to the partially concave surfaces.

Indications

- Finishing and polishing of aesthetic restorations with difficult to access proximal areas
- Aesthetic restorations with narrow interdental space conditions
- Excess (filling-cements, overhangs) removal on restoration surfaces

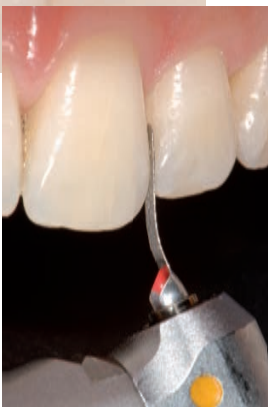

Clinical pictures by University of Zurich, Switzerland

Proxoshape Flexible						
L mm	8,5	8,5	8,5	11,0	11,0	11,0
µm	60	40	15	60	40	15
514	PSF60			PSF60L		
514		PSF40			PSF40L	
504			PSF15			PSF15L

To be used in combination with KaVo 61 LRG or KaVo contra-angle Intracompact 2061 CHC with 0.8 mm stroke length.